



Follicle aspiration needle

Knowledge for life since 1951.

Follicle aspiration needle



Knowledge for life since 1951.

Follicle aspiration needle

TIK d.o.o. | info@tik.si | +386 05 389 07 00


Knowledge for life since 1951. 31



Follicle aspiration needle

IVF -TIK

- Follicle aspiration needle is used for ultrasonic guided, transvaginal harvesting of oocytes from follicles
- The internal surface of the needle and tubing are extremely smooth which allows safe way for eggs to the test tube
- Enabled flushing of needle and connecting tube after the procedure
- Single lumen

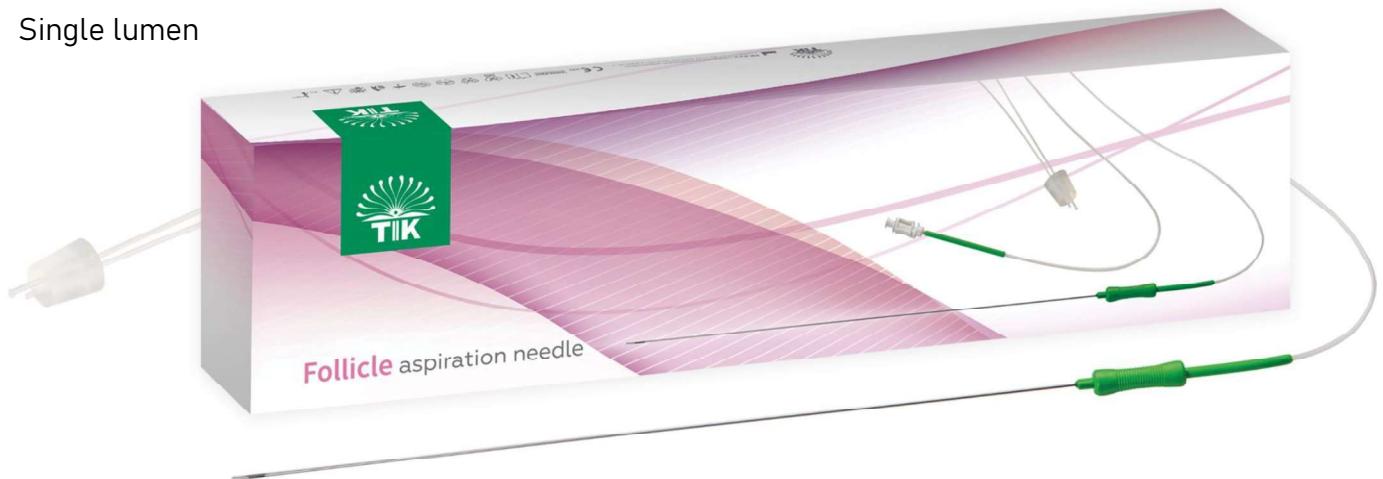

Follicle aspiration needle



Follicle aspiration needle



- Only biocompatible materials used
- Each needle is individually packed in a sterile packaging to avoid contamination up to the point of use
- To meet high standards, every needle before release is inspected by microscopical inspection to ensure compliance and at the same time provide consistent quality
- In order to provide Lot-to-lot consistency assembly and packaging activities are carried out within a certified quality management system (from raw material to finished product), at our premises





Follicle aspiration needle

EASY TO HANDLE:

- Ergonomically designed needle holder for easy control and orientation of the needle bevel. Needle holder is made of biocompatible material free of plasticizers(DEHP free)
- The thumb notch indicates the orientation of the bevel
- Additional sleeve prevents kinking of tubes

ECHOMARKING FOR OPTIMAL CONTROL:

- The echo marking enables high reflection and visibility for precise guidance with ultrasound used during oocyte retrieval.



Follicle aspiration needle



Follicle aspiration needle

CANNULA MADE FROM STAINLESS STEEL:

- Sharp needle bevel allows precision puncture using minimal force, high precision.

TWO BEVEL TYPES:

- Lancet and back cut bevel

SILICONE PLUG:

- Made from soft transparent silicone to secure good connection to the test tube



Follicle aspiration needle



Follicle aspiration needle

VACUUM TUBING AND ASPIRATION TUBING:

- Vacuum and aspiration pump tubing with different lengths
- Transparent, low friction tubing made of PTFE prevents oocyte damage between the needle and the test tube
- Funnel-shaped PTFE tubing endings enable safe outflow into collection tube

LUER LOCK CONNECTOR:

- Female luer lock connector enables connecting to an aspirator or a syringe
- Additional reinforcement of PTFE tubes at luer lock connector prevents kinking of tubes





Follicle aspiration needle

REF	OUTER DIAMETER OF THE NEEDLE (G/mm)	INTERNAL DIAMETER OF THE NEEDLE (mm)	NEEDLE LENGTH (mm)	LENGTH OF TUBING FROM NEEDLE TO PLUG (mm)	LENGTH OF TUBING FROM PLUG TO ASPIRATOR (mm)	ECHO MARKER	ADDITIONAL NEEDLE FOR FLUSHING	BEVEL TYPE		
58831	17/1,4	1,18	350	670	700	✓	✗	Lancet	10	120
59123	17/1,5	1,20	320	970	500	✓	✗	Lancet	10	120
59425	17/1,5	1,20	320	670	700	✗	✗	Lancet	10	120
59436	17/1,5	1,20	320	670	700	✗	✓	Lancet	10	120
59587	17/1,5	1,20	240	670	700	✓	✗	Lancet	10	120
59588	17/1,5	1,20	320	670	700	✓	✗	Lancet	10	120
63393	17/1,4	0,92	320	970	1000	✓	✗	Back cut bevel	10	120
63781	17/1,4	0,92	320	670	700	✓	✗	Back cut bevel	10	120

Follicle aspiration needle